# **INSTRUCTION MANUA**

# H197705

Silica Low Range Photometer





# Dear Customer,

Thank you for choosing a Hanna Instruments product.

Please read this instruction manual carefully before using this instrument. This manual will provide you with the necessary information for correct use of this instrument, as well as a precise idea of its versatility. If you need additional technical information, do not hesitate to e-mail us at tech@hannainst.com or view our worldwide contact list at www.hannainst.com.

# **TABLE OF CONTENTS**

1. PRELIMINARY EXAMINATION	4
2. SAFETY MEASURES	5
3. ABBREVIATIONS	5
4. SPECIFICATIONS	6
5. DESCRIPTION	7
5.1. GENERAL DESCRIPTION & INTENDED USE	7
5.2. FUNCTIONAL DESCRIPTION	8
5.3. PRECISION & ACCURACY	9
5.4. PRINCIPLE OF OPERATION	9
5.5. OPTICAL SYSTEM	10
6. GENERAL OPERATIONS	11
6.1. METER VALIDATION: CAL CHECK & CALIBRATION	11
6.2. CHEMICAL FORMULA & UNIT CONVERSION	12
6.3. GLP	13
6.4. LOGGING DATA & LOG RECALL	13
6.5. GENERAL SETUP	
6.6. REAGENTS & ACCESSORIES	16
6.7. CONTEXTUAL HELP	16
6.8. BATTERY MANAGEMENT	17
7. PHOTOMETER	18
7.1. COLLECTING & MEASURING REAGENTS AND SAMPLES	18
7.2. CUVETTE PREPARATION	18
8. METHOD PROCEDURE	20
9. WARNING & ERROR DESCRIPTIONS	23
10. BATTERY REPLACEMENT	25
11. ACCESSORIES	25
11.1. REAGENT SETS	25
11.2. OTHER ACCESSORIES	25
CERTIFICATION	26
RECOMMENDATIONS FOR USERS	26
WARRANTY	26

### 1. PRELIMINARY EXAMINATION

Remove the instrument and accessories from the packaging and examine it carefully. For further assistance, please contact your local Hanna Instruments Office or email us at tech@hannainst.com.

Each H197705C is delivered in a rugged carrying case and is supplied with:

- Sample cuvette (2 pcs.)
- Sample cuvette cap (2 pcs.)
- Plastic stopper (2 pcs.)
- A ZERO CAL Check Cuvette A
- H197705B CAL Check Cuvette B for Silica Low Range
- Cloth for wiping cuvettes
- Scissors
- 1.5V AA Alkaline batteries
- CAL Check standard certificate
- Instrument quality certificate
- Instruction manual

Each H197705 is delivered in a cardboard box and is supplied with:

- Sample cuvette (2 pcs.)
- Sample cuvette cap (2 pcs.)
- Plastic stopper (2 pcs.)
- 1.5V AA Alkaline batteries
- Instrument quality certificate
- Instruction manual

**Note:** Save all packing material until you are sure that the instrument works correctly. Any damaged or defective item must be returned in its original packing material with the supplied accessories.

### 2. SAFETY MEASURES



- The chemicals contained in the reagent kits may be hazardous if improperly handled.
- Read the Safety Data Sheets (SDS) before performing tests.
- Safety equipment: Wear suitable eye protection and clothing when required, and follow instructions carefully.
- Reagent spills: If a reagent spill occurs, wipe up immediately and rinse
  with plenty of water. If reagent contacts skin, rinse the affected area
  thoroughly with water. Avoid breathing released vapors.
- Waste disposal: For proper disposal of reagent kits and reacted samples, contact a licensed waste disposal provider.

# 3. ABBREVIATIONS

mg/L milligrams per liter (ppm)

mL milliliter

°C degree Celsius
°F degree Fahrenheit

GLP Good Laboratory Practice

HDPE High Density Polyethylene

LED Light Emitting Diode

LR Low Range

NIST National Institute of Standards and Technology

# 4. SPECIFICATIONS

Silica LR	Range	0.00 to 2.00 mg/L (as $SiO_2$ )
	Resolution	0.01 mg/L
	Accuracy	$\pm0.03$ mg/L $\pm3\%$ of reading at 25 °C
	Method	Adaptation of the ASTM Manual of Water and Environmental Technology, D859, Heteropoly Molybdenum Blue Method
	Light source	Light Emitting Diode
	Bandpass filter	610 nm
Measurement System	Bandpass filter bandwidth	8 nm
	Bandpass filter wavelength accuracy	±1.0 nm
	Light detector	Silicon photocell
	Cuvette type	Round 24.6 mm diameter (22 mm inside)
	Auto logging	50 readings
	Display	128 x 64 pixel B/W LCD with backlight
Additional Specifications	Auto-off	After 15 minutes of inactivity (30 minutes before a READ measurement)
	Battery type	1.5 V AA Alkaline (3 pcs.)
	Battery life	> 800 measurements (without backlight)
	Environment	0 to 50 °C (32 to 122 °F); 0 to 100% RH, non-serviceable
	Dimensions	142.5 x 102.5 x 50.5 mm (5.6 x 4.0 x 2.0")
	Weight (with batteries)	380 g (13.4 oz.)
	Case ingress protection rating	IP67, floating case

### 5. DESCRIPTION

### 5.1. GENERAL DESCRIPTION & INTENDED USE

The H197705 is an auto-diagnostic portable photometer that benefits from Hanna's years of experience as a manufacturer of analytical instruments. It has an advanced optical system that uses a Light Emitting Diode (LED) and a narrow band interference filter that allows for accurate and repeatable readings.

The optical system is sealed from outside dust, dirt and water. The meter uses an exclusive positive-locking system to ensure that the cuvettes are placed into the holder in the same position every time.

With the CAL Check  $^{\text{TM}}$  functionality, users are able to validate the performance of the instrument at any time and apply a user calibration (if necessary). Hanna Instruments CAL Check cuvettes are made with NIST traceable standards.


The built-in tutorial mode guides users step-by-step through the measurement process. It includes all steps required for sample preparation, the required reagents and quantities.


The H197705 meter measures the silica content in water samples in the 0.00 to 2.00 mg/L (ppm). The method is an adaptation of the ASTM Manual of Water and Environmental Technology, D859, Heteropoly Molybdenum Blue Method.

The H197705 photometer is a compact and versatile meter suitable for field or bench measurements, featuring a:

- Sophisticated optical system
- Meter validation using certified CAL Check cuvettes
- Tutorial mode guides the user step-by-step
- Auto logging
- Waterproof IP67, floating case
- GLP features

### 5.2. FUNCTIONAL DESCRIPTION

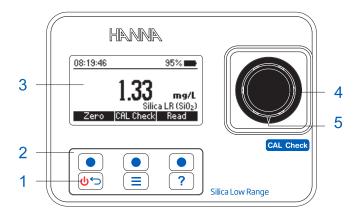

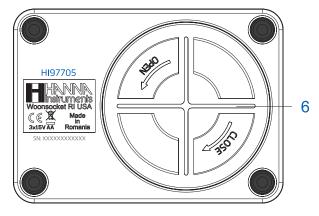

- 1. ON/OFF power button
- 3. Liquid Crystal Display (LCD)
- 5. Indexing mark

- 2. Keypad
- 4. Cuvette holder
- 6. Battery cover

# **Keypad Description**

The keypad contains 3 direct keys and 3 functional keys with the following functions:



Press the functional key to perform the function displayed above it on the LCD.



Press and hold to power off/on. Press briefly to return to the previous screen.



Press to access the menu screen.




Press to display the context-sensitive help menu.

### 5.3. PRECISION & ACCURACY

Precision is how closely repeated measurements are to one another. Precision is usually expressed as standard deviation (SD).

Accuracy is defined as the closeness of a test result to the true value.

Although good precision suggests good accuracy, precise results can be inaccurate. The figure explains these definitions.

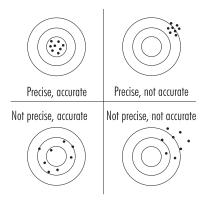

### 5.4. PRINCIPLE OF OPERATION

Absorption of light is a typical phenomenon of interaction between electromagnetic radiation and matter. When a light beam crosses a substance, some of the radiation may be absorbed by atoms, molecules or crystal lattices. Photometric chemical analysis is based on specific chemical reactions between a sample and reagent to produce a light-absorbing compound.

If pure absorption occurs, the fraction of light absorbed depends both on the optical path length through the matter and on the physical-chemical characteristics of the substance according to the Lambert-Beer Law. If all other factors are constant, the concentration "c" can be calculated form the absorbance of the substance.

Lambert-Beer Law:

-log 
$${\rm I/I_o} = \epsilon_\lambda\,{\rm c}\,{\rm d}$$
 or 
$${\rm A} = \epsilon_\lambda\,{\rm c}\,{\rm d}$$


 $I_o = intensity of incident light beam$ 

I = intensity of light beam after absorption $\epsilon_{\lambda} = \text{molar extinction coefficient at wavelength } \lambda$ 

c = molar concentration of the substance

d = optical path through the substance

### 5.5. OPTICAL SYSTEM

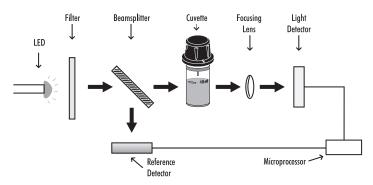

Instrument Block Diagram

The internal reference system (reference detector) of the H197705 photometer compensates for any drifts due to power fluctuations or ambient temperature changes, providing a stable source of light for your blank (zero) measurement and sample measurement.

LED light sources offer superior performance compared to tungsten lamps. LEDs have a much higher luminous efficiency, providing more light while using less power. They also produce little heat, which could otherwise affect electronic stability. LEDs are available in a wide array of wavelengths, whereas tungsten lamps have poor blue/violet light output.

Improved optical filters ensure greater wavelength accuracy and allow a brighter, stronger signal to be received. The end result is higher measurement stability and less wavelength error.

A focusing lens collects all of the light that exits the cuvette, eliminating errors from cuvette imperfections and scratches, eliminating the need to index the cuvette.

### 6. GENERAL OPERATIONS

### 6.1. METER VALIDATION: CAL CHECK & CALIBRATION


Validation of the H197705 involves verifying the concentration of the certified CAL Check standards. The CAL Check screen guides the user step-by-step through the validation process and user calibration (if necessary).

**WARNING:** Do not use any solutions or standards other than the Hanna Instruments  $^{\textcircled{o}}$  CAL Check Standards. For accurate validation and calibration results, please perform these at room temperature, 18 to 25  $^{\circ}$ C (64.5 to 77.0  $^{\circ}$ F).

**Note:** CAL Check Standards will not read the specified value in measurement mode. Protect the CAL Check cuvettes from direct sunlight by keeping them in the original packing. Store between 5 and 30  $^{\circ}$ C (41 to 86  $^{\circ}$ F), do not freeze.

To perform a CAL Check:

1. Press CAL Check from measurement mode.

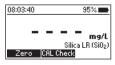

The "Not Available" message or the date, time and status of the last CAL Check will be displayed on the screen.





- 2. Press **Check** to start a new CAL Check. Press the bey at any time to abort the validation process.
- 3. Use the functional keys to enter the certificate value of the calibration standard found on the CAL Check Standard Certificate. Press **Next** to continue.



**Note:** This value will be saved in the instrument for future validation. If a new set of calibration standards is obtained, please update the certificate value.

4. Insert the A ZERO CAL Check Cuvette A then press **Next** to continue. The "Please wait..." message will be displayed during the measurement.





Insert the HI97705B CAL Check Cuvette B then press Next to continue. The "Please wait..." message will be displayed during the measurement.






- 6. When the CAL Check is complete the display will show one of the following messages and the value obtained during the measurement:
  - "PASSED": The measured value is within the accuracy specification, no user calibration is required.



"OUT OF SPECIFICATION" and Calibrate is available: The measured value
is near the expected value. To update the user calibration press Calibrate. Press
Accept to confirm or Cancel to return to the previous screen.



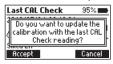



 "OUT OF SPECIFICATION": A user calibration is not allowed, the measured value is outside of the tolerance window. Check the certified value, expiration date and clean the outside of the cuvette. Repeat the CAL Check procedure. If this error continues, contact your nearest Hanna Instruments Customer Service Center.



### 6.2. CHEMICAL FORMULA & UNIT CONVERSION

Chemical formula and unit conversion factors are pre-programmed into the instrument. In order to view the displayed result in the desired chemical formula enter menu by pressing and use the functional keys to select *Chemical Form.* Press **Select** to change the displayed chemical formula. Use the functional keys to highlight the desired chemical formula and press **Select**. The selected formula will be saved when the instrument is powered off.





### 6.3. GLP

Press the \(\begin{align\*} \equiv \text{key to enter the menu.} Use the functional keys to select \(GLP\) and press \(Select\). Good Laboratory Practice (GLP) shows the date and time of the last user calibration (if available) or factory calibration. To erase the last user calibration and to clear the CAL Check press \(Clear\) and follow the prompts. Press \(Yes\) to erase and return to the factory calibration data or \(No\) to exit the clear procedure.



### 6.4. LOGGING DATA & LOG RECALL


The instrument features a data autolog function to help users keep track of all measurements. Every time a measurement is made the data is automatically saved. The data log can hold 50 individual measurements. When the data log is full (50 data points) the meter will rewrite the oldest data point.

Viewing and deleting the data is possible using the **Log Recall** menu.

Press the  $\equiv$  key to enter the menu. Use the functional keys to select *Log Recall* and press **Select**.




Use the functional keys to highlight a log and press **Info** to view additional information about the log. From this screen **Next** and **Previous** can be used to view other logs.

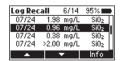



Press **Delete** to erase logged data. After pressing **Delete** a prompt on display is asking for confirmation.

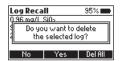



Press **No** or the by key to return to the previous screen. Press **Yes** to delete selected log.

Press **Del All** to erase all the logged data. If **Del All** is pressed, follow the prompt to confirm. Press **Yes** to delete all logged data, **No** or the beginning to the log recall.


### 6.5. GENERAL SETUP


Press the \( \equiv \) key to enter the menu. Use the functional keys to select *Setup* and press **Select**. Use the functional keys to highlight desired option.

# **Backlight**

Option: 0 to 100 %

Press **Modify** to access the backlight intensity. Use the functional keys to increase or decrease the value. Press **Accept** to confirm or the beyone to the setup menu without saving the new value.

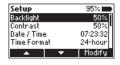


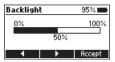


### Contrast

Option: 0 to 100 %

Press **Modify** to change the display's contrast. Use the functional keys to increase or decrease the value. Press **Accept** to confirm the value or the key to return to the *Setup* menu without saving the new value.

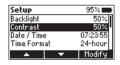


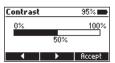

### Date & Time

Press Modify to change the date and time. Press the functional keys to highlight the value to be modified (year, month, day, hour, minute or second). Press Edit to modify the highlighted value. Use the functional keys to change the value.

Press **Accept** to confirm or the **b** key to return to the previous screen.

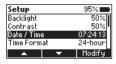





### Time Format


# Option: AM/PM or 24-hour


Press the functional key to select the desired time format.



### **Date Format**

Press **Modify** to change the date format. Use the functional keys to select the desired format. Press **Accept** to confirm or the **bestup** key to return to the **Setup** menu without saving the new format.

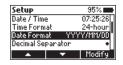

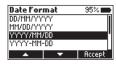

# **Decimal Separator**

Option: Comma (,) or Period (.)

Press the functional key to select the desired decimal separator. The decimal separator is used on the measurement screen.



# Language

Press **Modify** to change the language. Use the functional keys to select the desired language. Press **Accept** to choose one of the languages installed.






# Beeper

# Option: Enable or Disable


When enabled, a short beep is heard every time a key is pressed. A long beep alert sounds when the pressed key is not active or an error is detected. Press the functional key to enable or disable the beeper.

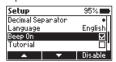

### Tutorial

### Option: Enable or Disable

When enabled, the user will be guided step-by-step through the measurement procedure.

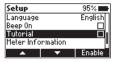

### Meter Information

Press **Select** to view the model, serial number, firmware version and selected language. Press the beginning to return to the *Setup* menu.





### **Restore Factory Settings**

Press **Select** to reset to factory settings. Press **Accept** to confirm or **Cancel** to exit without restoring the factory settings.





### 6.6. REAGENTS & ACCESSORIES

Press the (=) key to enter the menu. Use the functional keys to select *Reagents / Accessories* and press **Select** to access a list of reagents and accessories. To exit press the (0) key.






### 6.7. CONTEXTUAL HELP

The H197705 offers an interactive contextual help mode that assists the user at any time. To access the help screen press the ? key.



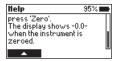


The instrument will display additional information related to the current screen. To read all the available information, scroll the text using the functional keys.

To exit help mode press the  $\bigcirc$  or the  $\bigcirc$  key and the meter will return to the previous screen.


### 6.8. BATTERY MANAGEMENT

The meter will perform an auto-diagnostic test when it is powered on. During this test, the Hanna Instruments  $^{\circledR}$  logo will appear on the LCD. If the auto-diagnostic test was successful, the meter is ready for use. The battery icon on the LCD will indicate the battery status:


• battery full

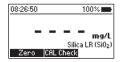

• battery below 10%, replace the batteries soon

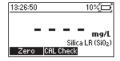

• battery is low, replace the batteries with new ones

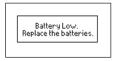

To conserve battery, the meter will turn off automatically after 15 minutes of inactivity. If a zero reading has been done but not a read, auto-off time is increased to 30 minutes.

### 7. PHOTOMFTER

### 7.1. COLLECTING & MEASURING REAGENTS AND SAMPLES


# **Proper Use of Powder Packet**

- 1. Use scissors to open the powder packet.
- 2. Push the edges of the packet to form a spout.
- 3. Pour out the content of the packet.



# **Proper Use of Dropper Bottle**

- 1. Tap the dropper on the table several times and wipe the outside of the tip with a cloth.
- 2. Always keep the dropper bottle in a vertical position while dosing the reagent.

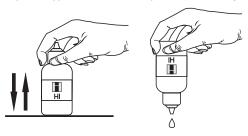

### 7.2. CUVETTE PREPARATION

Proper mixing is very important for reproducibility of the measurements. The proper mixing technique is listed in the method procedure.

(a) The mixing method is indicated with "shake gently" using one of the following icons:



(b) The mixing method is indicated with "swirl" using one of the following icons:



In order to avoid reagent leaking and to obtain more accurate measurements, close the cuvette first with the supplied HDPE plastic stopper and then the black cap.



Whenever the cuvette is placed into the measurement holder, it must be dry outside and free of fingerprints, oil or dirt. Wipe it thoroughly with HI731318 microfiber cleaning cloth or a lint-free wipe prior to insertion.



Shaking the cuvette can generate bubbles in the sample, causing higher readings. To obtain accurate measurements, remove such bubbles by swirling or by gently tapping the cuvette.

Do not let the reacted sample stand too long after reagent has been added. For best accuracy, respect the timings described in the method.

It is possible to take multiple readings in a row, but it is recommended to take a new zero reading for each sample and to use the same cuvette for zeroing and measurement when possible.

Discard the sample immediately after the reading has been taken, or the glass might become permanently stained.

All the reaction times reported in this manual are at 25 °C (77 °F). In general, the reaction time should be increased for temperatures lower than 20 °C (68 °F), and decreased for temperatures higher than 25 °C (77 °F).

# 8. METHOD PROCEDURE

### REQUIRED REAGENTS

Code	Description	Quantity
H193705A-0	Silica LR Reagent A	6 drops
H193705B-0	Silica LR Reagent B	1 packet
H193705C-0	Silica LR Reagent C	1 packet

### **REAGENT SETS**


HI93705-01	Silica LR Reagent - 100 tests
HI93705-03	Silica LR Reagent - 300 tests
For other accesso	ries see ACCESSORIES section.

### MEASUREMENT PROCEDURE


**Note:** If tutorial mode is disabled, follow the measurement procedure below. If the tutorial mode is enabled, press **Measure** and follow the messages on the screen.

• Fill the cuvette with 10 mL of unreacted sample (up to the mark).

 Add 6 drops of HI93705A-0 Silica LR Reagent A. Replace the plastic stopper and the cap. Swirl to mix the solution.

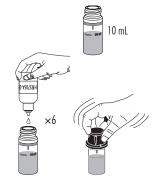

 Press Zero. The display will show a 4 minute countdown prior to adding H193705B-0 Silica LR Reagent B. To skip the timer, press Zero. To pause the 4 minute timer, press Stop. To start the 1 minute timer, press 1:00 min.

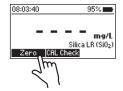



Add the content of one packet of HI93705B-0
 Silica LR Reagent B. Replace the plastic
 stopper and the cap. Shake gently until it is
 completely dissolved.



 Press Start. The display will show a 1 minute countdown prior to zero. To skip the timer, press Zero.




• Insert the cuvette into the holder and ensure that the notch on the cap is positioned securely in the groove.




 Press Zero. The display will show "-0.0-" when the meter is zeroed and ready for measurement.



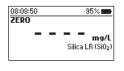




- Remove the cuvette.
- Add the content of one packet of HI93705C-0 Silica LR Reagent C. Replace the plastic stopper and the cap. Shake gently until it is completely dissolved.

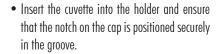






Press Read. The display will show a 3 minute countdown prior to the measurement.
 To skip the timer, press Read twice. When the timer ends, the meter will perform the reading. The instrument displays the results in mg/L of silica (SiO<sub>2</sub>).

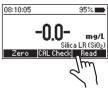





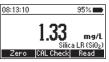

• Press the  $(\equiv)$  key and use the functional keys to select *Chemical Form*.





 Use the functional keys and press Select to change the displayed chemical formula to mg/L of silicon (Si).





### **INTERFERENCES**

Interference may be caused by:

- $\bullet$  Phosphate above 75 mg/L, causes an 11% reduction in reading
- Phosphate above 60 mg/L, causes a 2% reduction in reading
- Sulfide and high concentration of iron
- Eliminate color and turbidity interferences by zeroing the meter with the original water sample

# 9. WARNING & ERROR DESCRIPTIONS

The instrument shows clear warning messages when erroneous conditions appear and when measured values are outside the expected range.

The information below provides an explanation of the errors and warnings, and recommended action to be taken.



There is an excess amount of ambient light reaching the detector. Ensure that the notch on the cap is positioned securely in the groove before performing any measurements. If the issue persists, please contact Hanna Instruments technical support.



The sample and the zero cuvettes are inverted. Swap the cuvettes and repeat the measurement.



There is either too much light or the instrument can not adjust the light level. Please check the preparation of the zero cuvette and that the sample does not contain any debris.



under limit. Put the meter in a warm place. Continue

Warning

The meter is either overheating or its temperature has dropped too low to operate within published accuracy specifications. The meter must be between 0 and 50  $^{\circ}$ C (32 and 122  $^{\circ}$ F) to perform any measurements.



Meter temperature has changed significantly since the zero measurement has been performed. The zero measurement must be performed again.




The measured value is outside the limits of the method. Verify that the sample does not contain any debris. Check the sample preparation and the measurement preparation.




Date and time settings have been lost. Please reset the values. If the issue persists, please contact Hanna Instruments technical support.



English is the only available language. Help function is not available. Restart the meter. If the issue persists, please contact Hanna Instruments technical support.

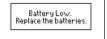

Battery level is too low for the meter to function properly. Replace the batteries with new ones.

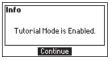

Tutorial mode has been enabled in the Setup menu. Press **Continue** and follow the prompt on the screen. Tutorial mode can be disabled in the Setup menu.




A critical error has occurred. Restart the meter. If the issue persists, please contact Hanna Instruments technical support.

# 10. BATTERY REPLACEMENT

To replace the instrument's batteries, follow these steps:

- Turn the instrument off by pressing and holding the 😊 key.
- Remove the battery cover by turning it counterclockwise.
- Remove the old batteries, replace them with three new 1.5V AA batteries.
- Replace the battery cover, turn it clockwise to close.

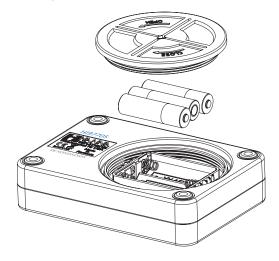

# 11. ACCESSORIES

# 11.1. REAGENT SETS

Code	Description
HI93705-01	Silica LR Reagent - 100 tests
HI93705-03	Silica LR Reagent - 300 tests

# 11.2. OTHER ACCESSORIES

Code	Description
HI7101412	blue carrying case for H1977xx and 2 CAL Check cuvettes
HI731318	cloth for wiping cuvettes (4 pcs.)
HI731331	glass cuvette (4 pcs.)
HI731336N	cap for glass cuvette (4 pcs.)
HI93703-50	cuvette cleaning solution (230 mL)
HI97705-11	CAL Check <sup>®</sup> standards for Silica LR - cuvette kit

### **CERTIFICATION**

All Hanna Instruments conform to the CE European Directives.



**Disposal of Electrical & Electronic Equipment**. The product should not be treated as household waste. Instead hand it over to the appropriate collection point for the recycling of electrical and electronic equipment which will conserve natural resources. **Disposal of waste batteries**. This product contains batteries, do not dispose of them with other household waste. Hand them over to the appropriate collection point for recyclina.

Ensuring proper product and battery disposal prevents potential negative consequences for the environment and human health. For more information, contact your city, your local household waste disposal service, the place of purchase or go to www.hannainst.com.



# RECOMMENDATIONS FOR USERS

Before using this product, make sure it is entirely suitable for your specific application and for the environment in which it is used. Any variation introduced by the user to the supplied equipment may degrade the meter's performance. For yours and the meter's safety do not use or store the meter in hazardous environments.

# WARRANTY

The H197705 is warranted for two years against defects in workmanship and materials when used for its intended purpose and maintained according to instructions. This warranty is limited to repair or replacement free of charge. Damage due to accidents, misuse, tampering or lack of prescribed maintenance is not covered.

If service is required, contact your local Hanna Instruments Office. If under warranty, report the model number, date of purchase, serial number (engraved on the bottom of the meter) and the nature of the problem. If the repair is not covered by the warranty, you will be notified of the charges incurred. If the meter is to be returned to Hanna Instruments, first obtain a Returned Goods Authorization (RGA) number from the Technical Service department and then send it with shipping costs prepaid. When shipping any meter, make sure it is properly packed for complete protection.



# World Headquarters

Hanna Instruments Inc. Highland Industrial Park 584 Park East Drive Woonsocket, RI 02895 USA www.hannainst.com



MAN97705 Printed in ROMANIA